

Technical data sheet

# K-365 Dist Line



## Dist Line

Due to its modular concept the Dist systems can be adapted to various needs. The different instrument configurations allow to meet different requirements.

### A) Three instrument performance levels

- EasyDist: Entry level unit
- BasicDist: Intermediate unit
- MultiDist: Advanced unit

### B) Splash protectors

Advantage of the glass splash protector:

- visual monitoring of the distillation process makes method development more convenient

Advantage of the plastic splash protector:

- increased durability to sodium hydroxide and extended lifetime

### C) Application versatility

All Dist system are designed to distill steam volatile compounds such as:

- Alcohol
- Essential oils
- VDK (Vicinal diketones)

All Basic and Multi systems contain an acid resistant pump, which allows to distill acid released volatiles:

- Sulfite
- Cyanide
- Volatile acids
- Phenol
- Formaldehyde

The Multi systems additionally contain a sodium hydroxide resistant pump for Kjeldahl applications:

- Ammonia
- Total
- Kjeldahl Nitrogen (TKN)
- Protein Nitrate/Nitrite

### D) External titrators

Various titrators can be connected to the MultiDist systems:

- Metrohm (Eco Titrator and TitrinoPlus 848/877)
- SI Analytics (TitroLine 5000, TitroLine 7000 and Titronic 300)
- Mettler Toledo (G10S, G20S, T5, T7 and T9)

All other titrator models can be used via a semi-automated integration.

## Scope of delivery

All configurations are supplied ready to use.

	EasyDist	BasicDist	MultiDist
Instrument	1	1	1
Splash protector	1*	1*	1*
Titration set	-	-	*
External titrator	-	-	*
Recirculating chiller	*	*	*
Power cable	1	1	1
Sample tube (300 mL)	1	1	1
Tanks (10 L)	1	2	3
Cap Set, 10 L tank	1	2	3
Tanks (20 L)	-	-	1
Cap Set, 20 L tank	-	-	1
Suction tube tank	1	2	5
Hose for chemical feed (2 m)	1	1	4
Hose for cooling water feed (1.5 m)	1	1	1
Hose for cooling water drain (3 m)	1	1	1
Pair of glass tongs	1	1	1
Hose for instrument drain	-	*	0
Hose clip acid pump	-	1	1
Hose clip reagent pumps	-	1	4
Hose clip peristaltic pumps	-	*	2

\* According to order code

## Order code EasyDist

Choose the configuration according to your needs:

**1 1 K 3 6 5 1 2 0**    

Splash Protector

- 4 Glass without Dosing Inlet
- 5 Plastic without Dosing Inlet

Recirculation Chiller

- 0 None
- 1 F-314 230V

## Order code BasicDist

Choose the configuration according to your needs:

**1 1 K 3 6 5 2 2**      

Configuration

- 0 Base Unit
- 1 Base Unit with Sample Tube Aspiration
- 2 Base Unit with Sample Tube Aspiration and I-300 Pro

Splash Protector

- 1 Glass with Dosing Inlet
- 2 Plastic with Dosing Inlet

Recirculation Chiller

- 0 None
- 1 F-314 230V

## Order code MultiDist

Choose the configuration according to your needs:



### Configuration

- 0 Base Unit
- 1 Base Unit with Titration Setup and Receiver Waste Aspiration
- 2 Base Unit with Titration Setup, Receiver Waste Aspiration and Metrohm Eco Titrator (potentiometric)

### Splash Protector

- 1 Glass with Dosing Inlet
- 2 Plastic with Dosing Inlet
- 3 Devarda

### Recirculation Chiller

- 0 None
- 1 F-314 230V

## Technical data

### K-365 Dist Line

	EasyDist	BasicDist	MultiDist
Dimensions (W x D x H)	320 x 400 x 730 mm	320 x 400 x 730 mm	320 x 400 x 730 mm
Weight	23 kg	23 kg	23 kg
Connection voltage	220 - 240 ± 10 % VAC	220 - 240 ± 10 % VAC	220 - 240 ± 10 % VAC
Power consumption	2100 W	2100 W	2100 W
Frequency	50 / 60 Hz	50 / 60 Hz	50 / 60 Hz
IP Code	IP20	IP20	IP20
Pollution degree	2	2	2
Overvoltage category	II	II	II
Tap water pressure	1 - 10 bar	1 - 10 bar	1 - 10 bar
Tap water condition <small>(If the temperature is more than 25°C use a recirculating chiller)</small>	1- 25°C	1 - 25°C	1 - 25°C
Tap water consumption	~ 1.2 L/min	~ 1.2 L/min	~ 1.2 L/min
Approvals	CE, CSA	CE, CSA	CE, CSA

## Analysis specifications

	EasyDist	BasicDist	MultiDist
Measuring range	analyte specific	analyte specific	0.02 – 220 mg N
Analysis time for 30 mg N <small>(with Metrohm Eco Titrator without titrator preparation)</small>	-	-	3.5 min
Analysis time for 200 mg N <small>(with Metrohm Eco Titrator without titrator preparation)</small>	-	-	5 min
LOD (Limit of Detection) of N	-	-	0.007 mg
LOQ (Limit of Quantification) of N	-	-	0.02 mg
LOD (Limit of Detection) of SO <sub>2</sub>	-	-	0.026 mg
LOQ (Limit of Quantification) of SO <sub>2</sub>	-	-	0.078 mg
Recovery	analyte specific	analyte specific	N: >98% with digestion process N: >99.5% (direct distillation)
Measurement reproducibility <small>(Direct distillation of 1 mg N absolute and 300 sec. distillation )</small>	analyte specific	analyte specific	N: <0.8%

## Instrument specifications

	EasyDist	BasicDist	MultiDist
Distillation capacity (100%)	~ 40 mL/min	~ 40 mL/min	~ 40 mL/min
Distillation capacity (10%)	~ 12.5 mL/min	~ 12.5 mL/min	~ 12.5 mL/min
Adjustable distillation power	10 - 100%	10 - 100%	10 - 100%
MaxAccuracy Mode <small>(Compensation of power grid fluctuations by voltage stabilizer)</small>	Yes	Yes	Yes
Distillate reproducibility (RSD) <small>(At 300 seconds distillation time)</small>	<1%	<1%	<1%
AutoDist Mode <small>(Automatic detection of the condensation start)</small>	No	Yes	Yes
NaOH pump	No	No	Yes
H <sub>2</sub> O pump	No	No	Yes
Acid pump	No	Yes	Yes
Boric acid pump	No	No	Yes
Automatic sample disposal	No	Optional	Yes
Automatic receiver disposal	No	No	Optional
Safety sensors (safety door sensor, cooling water temperature sensor, sample tube in-place sensor, steam generator overpressure sensor)	Yes	Yes	Yes
Titration vessel	No	No	Optional
Titrator connection	No	No	Optional
IQ/OQ	No	Yes	Yes

## Interface specifications

	Interface EasyDist, BasicDist	Interface pro BasicDist (option), MultiDist
Display	LCD color, 4.3 in	LCD color, 7.0 in
Control elements	Buttons / control knob	Touch-screen / control knob
Max. number of methods	8	96
Max. number of determinations	36	512
Max. number of determinations per series	16	32
Max. number of templates	4	16
Max. number of users	4	32
Max. number of determination results	40	256
Language	en, de, fr, it, es, pt, ja, zh, ru, pl, ko, id	en, de, fr, it, es, pt, ja, zh, ru, pl, ko, id

## Connections

	EasyDist	BasicDist	MultiDist
LAN port	-	-	1
BUCHI standard communication ports (COM)	2	2	3
Stirrer port	-	-	1
Level detection sensor port	1	1	1
H <sub>2</sub> O tank level sensor port	1	1	1
H <sub>3</sub> BO <sub>3</sub> tank level sensor port	-	-	1
Acid tank level sensor port	-	-	-
NaOH tank level sensor port	1	1	1
Waste tank level sensor port	-	option	2
Dispenser port (RS-232)	-	-	1
USB port	-	option	1
USB port	-	option	1

## Ambient conditions

For indoor use only.

Max. altitude above sea level	2000 m
Ambient and storage temperature	5–40 °C
Maximum relative humidity	80% for temperatures up to 31 °C decreasing linearly to 50 % relative humidity at 40 °C

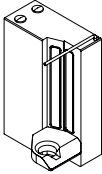
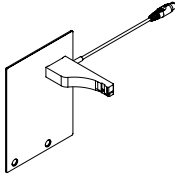
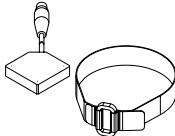
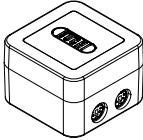
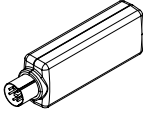

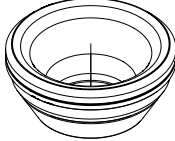
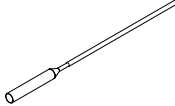
## Materials

Component	Materials of construction
Housing	Polyurethane
Housing electric	Stainless steel
Steam generator housing	Stainless steel
Housing stirrer	Glass fiber reinforced polypropylene
Rear cover	Polyvinyl chloride
Glass parts	Borosilicate 3.3
Plastic splash protector	Glass fiber reinforced polypropylene
Protective cover	PMMA
Sealing	CSM
NaOH hose	EPDM
Three-way valve	PP / PE

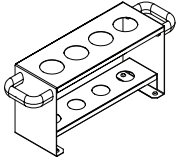
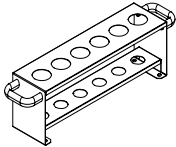
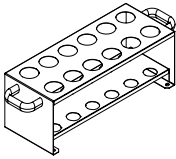

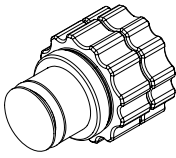


## Spare parts and accessories


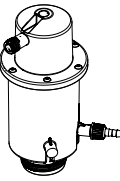
### Accessories

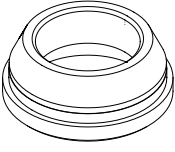
	Order no.	Image
Reaction detection sensor Detects the optimal alkalization point, saves up to 30 % of NaOH in every determination, saves resources and reduces the ecological footprint	11072666	
OnLevel sensor Sensor to set the desired distilled volume in the receiving vessel, for most reproducible results in every determination	11070270	
Tank level sensor kit Detects the current reagent / water / waste level in the tank, for higher convenience and smooth operation during runs	11072294	
Eco Titrator Including pre-defined methods for MultiKjel / MultiDist and connection to enable automatic titration	11072748	
StatusLight cpl., incl. communication cable Indicates the status of the instrument (instrument is ready to use, has an error or is in operation).	11068959	
BUCHI Bluetooth® Dongle, connects instrument to smartphone via Bluetooth®	11067770	
Adapter for 3rd party sample tubes Adapter setup for other sample tubes (250 mL), sealing included	11072398	
Sealing for sample tube 3rd party adapter	11072180	
Condenser outlet for alcohol application	043096	

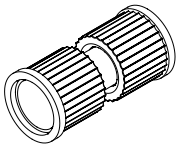
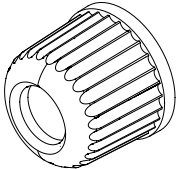

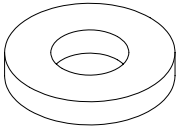
## Racks and rack parts

	Order no.	Image
Sample tube holder (4x 500 mL tubes)	016951	
Sample tube holder (6 x 300 mL tubes)	043039	
Sample tube holder (12 x 300 mL tubes)	043041	
SO <sub>2</sub> absorption glass set	11073599	
Cyanide caps Caps for hermetical closing of all types of Kjeldahl sample tubes during sample preparation for cyanide in food	11067871	
Dispenser unit Titronic 300 For back titration, including connection cable to distillation unit	11075111	

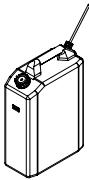
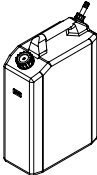
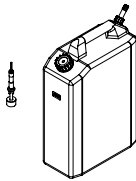
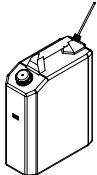
## Splash protector spare parts

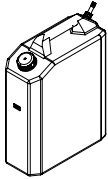

	Order no.	Image
Glass splash protector (EasyDist)	11071015	
Plastic splash protector (EasyDist)	11072055	

	Order no.	Image
Glass splash protector For visible distillation processes.	11071013	
Plastic splash protector The most durable splash protector.	11070670	
Devarda splash protector Prevents from foaming into the condenser.	11071014	
Sealing sample tube to splash protector	11073674	
Sealing sample tube to splash protector (acid resistant)	11073675	
Steam inlet tubing to sample tube	043424	
Steam inlet tubing to sample tube (750 mL)	043119	
Connecting piece (CSM) Rubber piece between splash protector and steam outlet tube, al- kali-resistant.	019002	
Connecting piece EPDM Rubber piece between splash protector and steam outlet tube, acid-resistant.	11062478	

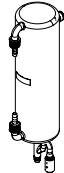
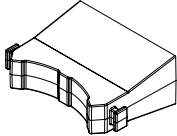
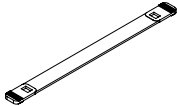
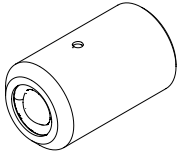
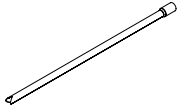
	Order no.	Image
Bridge splash protector to condenser	11070620	
Cap nut GL14 Inlets for splash protector	033577	
Hose barbs, set. 4 pcs, straight, GL 14, FPM seal Content: Hose barbs, cap nuts, seals.	040296	
Seals, set. 10 pcs, for hose barbs GL 14, FPM, green	040040	

## Tank spare parts


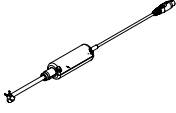
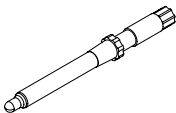
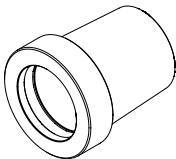
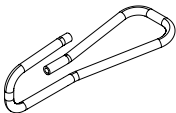
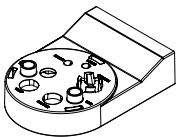
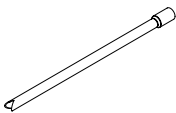
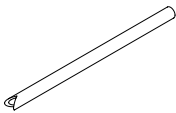
	Order no.	Image
Cap set, 10L tank	11073639	
Storage tank (10 L) for chemicals	043468	
Storage tank (10 L) for waste	043470	
Storage tank (10 L) with level sensor for waste	043474	
Cap Set, 20 L tank	11072174	
Storage tank (20 L) for chemicals Caps included	043469	

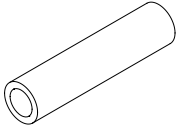
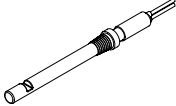
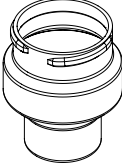
	Order no.	Image
Storage tank (20 L) for waste Caps included	043471	
Tank labels	043434	

## Condenser spare parts





	Order no.	Image
Glass condenser with check valve	11072183	
Clamp for condenser	11066868	
Rubber band Holding the glass condenser	11070669	
Check valve For overpressure safety	11071740	
Dist. Outlet PTFE L=300 for receiver side	11071940	

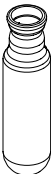


## Titration related parts

	Order no.	Image
Receiving vessel 420 mL	043390	
Stirrer For the receiving vessel.	11070246	
pH Electrode glass, refillable For potentiometric titrations.	11065834	
pH Electrode gel, maintenance-free For potentiometric titrations.	11056842	
Clamp cone pH-electrode	11069793	
Hose Tygon Ø 8.0 x 4.8 Hose to receiving module (waste aspiration, boric acid dosing)	043364	
Receiving support For attaching the receiving vessel to the distillation unit	11071003	
Condensate outlet to receiving module, PTFE	11071941	
Aspiration tubing receiving vessel	11072589	

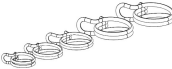
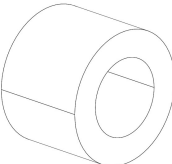
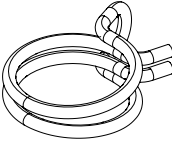


	Order no.	Image
H <sub>3</sub> BO <sub>3</sub> tubing receiving vessel For dosing H <sub>3</sub> BO <sub>3</sub>	11072637	
Colorimetric sensor (Optrode) For colorimetric titrations.	11066601	
Kit colorimetric titration (Optrode) Includes all Optrode accessories (w/o Optrode sensor and receiving vessel).	11068260	
Receiving vessel (Optrode) For colorimetric titrations with Optrode sensor.	11068263	

## Sample tubes

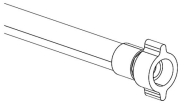
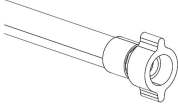
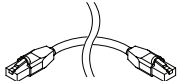
	Order no.	Image
Sample tubes micro (100 mL) 4 pcs.	11057442	
Sample tubes 300 mL 20 pcs.	11059690	
Sample tubes 300 mL For sample volumes up to 200 mL or 5 g in weight 4 pcs.	037377	
Sample tubes graduated 300 mL 4 pcs.	043049	
Sample tube 500 mL	026128	

	Order no.	Image
Sample tubes 500 mL 4 pcs.	043982	
Set sample tube 750 mL with long suction tubing	11058999	
Pair of glass tongs	002004	

## Cable and tubing

Specification	Order no.	Image
Hose peristaltic pump out (2.5 m)	11071630	
Hose clips peristaltic pump hose	043586	
Hose acid pump (2 m)	043588	
Hose clip D7.6	043571	
Hose to NaOH/ H <sub>2</sub> O/H <sub>3</sub> BO <sub>3</sub> pumps or drain (2 m)	11072687	
Hose clip to NaOH/H <sub>2</sub> O/ H <sub>3</sub> BO <sub>3</sub> pumps or drain	043841	
Suction to tube tank (0.58 m)	043407	



Specification	Order no.	Image
Tap water hose cpl.	037780	
Set of sealings tap water hose	040043	
Connection cable to SI-Analytics TitroLine/Titronic for K-365	043618	
Cable to Metrohm & Mettler (T-series) titrator	043617	
To connect the Metrohm Titrimo plus, a Metrohm remote box is needed, see also quick guide.		
Tubing to peristaltic pump	11070015	
Communication cable RJ45, 2 m	044989	
Connection between vacuum controller/interface and recirculating chiller, vacuum controller/interface and vacuum pump or steam distillation unit and Eco titrator.		

## Consumables

	Order no.
2% boric acid pH 4.65 +/- 0.15 with Sher indicator, 5 L	11064972
4% boric acid pH 4.65 +/- 0.15 with Sher indicator, 5 L	11064973
4% boric acid pH 4.65 +/- 0.15 with bromocresol green / methyl red mixed indicator, 5 L	11064976
Buffer solution pH 4.00 ± 0.01 twin-neck bottle, 250 mL	11064974
Buffer solution pH 7 ± 0.01 twin-neck bottle, 250 mL	11064975
KCl electrolyte, saturated	11064978
Electrode storage solution, 250 mL	
Sher indicator 100 mL Color indicator for colorimetric titrations.	003512
Ammonium phosphate monobasic, 25 g	045527
Reference substance for the verification and OQ of the Kjeldahl distillation units.	
Weighing boats	11060522
Nitrogen-free weighing paper 58 x 10 x 10 mm, 100 pcs.	
Kjeldahl Tablets Titanium	11072627
For efficient digestions, 1000 tablets (3.5 g K <sub>2</sub> SO <sub>4</sub> / 0.105 g CuSO <sub>4</sub> * 5 H <sub>2</sub> O / 0.105 g TiO <sub>2</sub> )	

Kjeldahl Tablets Titanium Micro 11072628

For small sample sizes, 1000 tablets (1.5 g  $K_2SO_4$  / 0.045 g  $CuSO_4 \cdot 5 H_2O$  / 0.045 g  $TiO_2$ )

Kjeldahl Tablets Missouri 11072629

For efficient digestions without Titanium, 1000 tablets (4.98 g  $K_2SO_4$  / 0.02 g  $CuSO_4 \cdot 5 H_2O$ )

Kjeldahl Tablets ECO 11072630

For eco-friendly digestions with minimum amount of copper, 1000 tablets (3.998 g  $K_2SO_4$  / 0.002  $CuSO_4 \cdot 5 H_2O$ )

Kjeldahl Tablets Antifoam 11072631

For foaming samples in digestions to avoid cross contamination due to foaming, 1000 tablets (0.97 g  $Na_2SO_4$  / 0.03 g Silicone)

Kjeldahl Tablets Copper Micro 11072632

For small sample sizes without Titanium, 1000 tablets (1.5 g  $K_2SO_4$  / 0.15 g  $CuSO_4 \cdot 5 H_2O$ )

## Maintenance kit

	Order no.
Customer Kit for Dist Line	11073025
Including suction tube with EPDM connecting piece, gaskets and EPDM seal with ring	

## Upgrade kits

We offer the following upgrade kits:

- Upgrade Basic base unit with aspiration pump
- Upgrade Basic base unit and aspiration pump with I-300 Pro
- Upgrade Basic base unit with aspiration pump and I-300 Pro
- Upgrade Multi base unit with titration vessel

For ordering these upgrades please contact your local service support.

## Documentation

	Order no.
IQ/OQ set MultiDist en Including Eco Titrator if available	11073602
Repeating OQ MultiDist en Including Eco Titrator if available	11073603
IQ/OQ set BasicDist en	11073631
Repeating OQ BasicDist en	11073632
Kjeldahl Knowledge Base (EN) Comprehensive guide covering all aspects of the theoretical and practical know-how.	11595478
Kjeldahl Practice Guide (EN) Provides theoretical background information, useful hints and calculation tables for daily routine work	11592548
Kjeldahl Practice Guide (DE)	11592547
Kjeldahl Practice Guide (CN)	11592549



11594202F en K-365 Dist Line 07.2024  
Technical data are subject to change without notice  
Quality Systems ISO 9001

---

We are represented by more than 100 distribution partners worldwide.  
Find your local representative at:

[www.buchi.com](http://www.buchi.com)

Quality in your hands

---