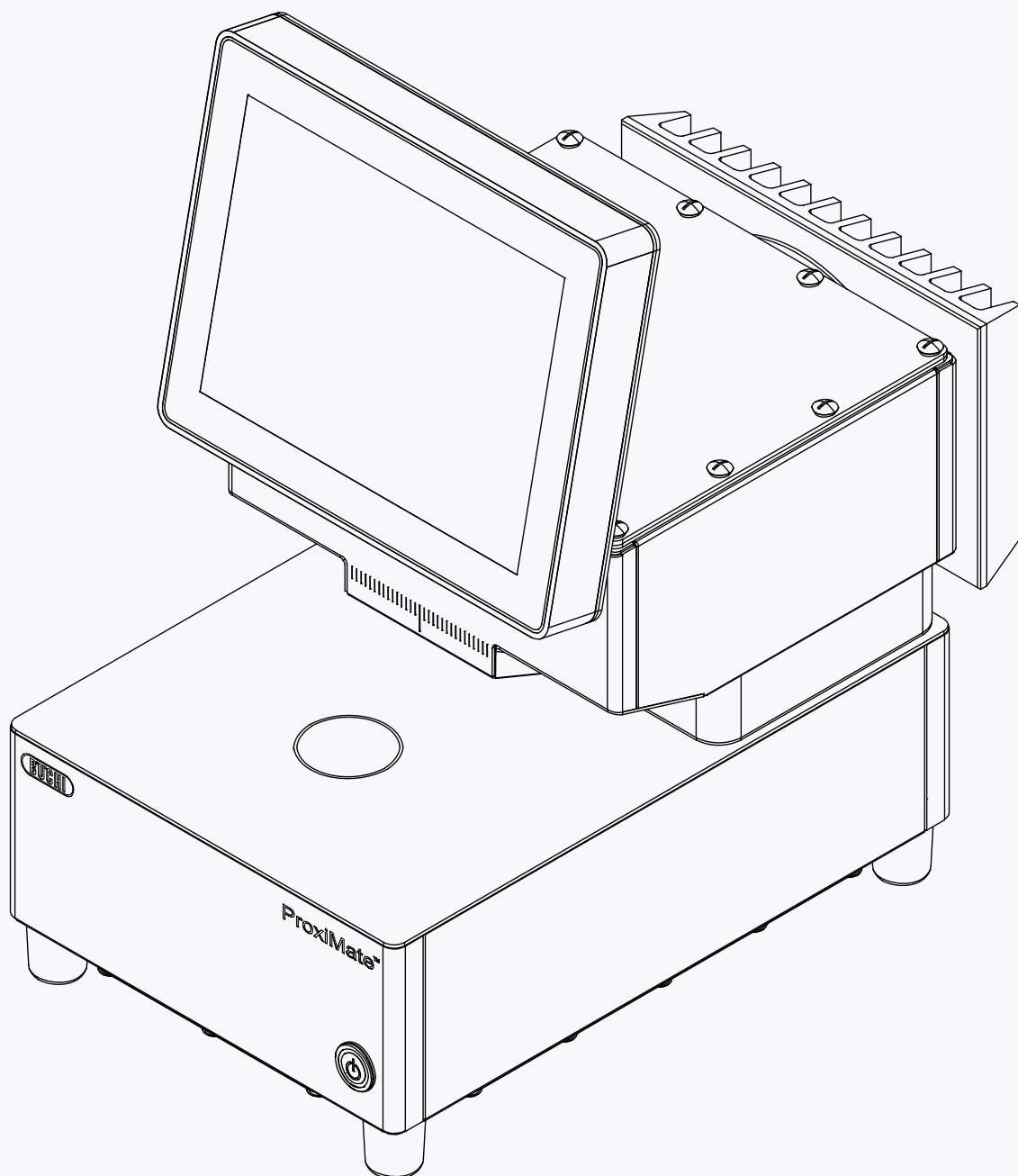




Technical data sheet

ProxiMate™

Depending on specification the BUCHI ProxiMate™ Essential is an NIR or visible light plus NIR spectrometer based on diode array technology. ProxiMate™ Essential is designed especially to work in the food and feed environments. The stainless steel case is designed to prevent corrosion when acidic, basic or chlorine detergents are applied.



Order code

Choose the configuration according to your needs:

	1	1	P	M						0	0	0	0
--	---	---	---	---	--	--	--	--	--	---	---	---	---

Spectrometer Type	
1 NIR	
2 VIS & NIR	

Optic Setup	
1 Up view	
2 Down view	
3 Dual view	

IP Module	
1 Standard	
2 Advanced	

Software Package	
1 NIRWise™ Operator Software	
2 NIRWise™ Operator Software including NIRWise™ PLUS Chemometrics Suite	

Accessories Sample Carrier	
1 Large	
2 Small *	

* Available with down view configuration only

Scope of delivery

All configurations are supplied ready to use.



NOTE

The scope of delivery depends on the configuration of the purchase order.

Accessories are delivered as per the purchase order, order confirmation, and delivery note.

Component	Specifications supplied
ProxiMate™	1
including white reference standard, spacers and positioning plate	
Operation Manual	1
Guide for electrical installation	1
Quick Guide	1

Technical data

ProxiMate™

Specification	ProxiMate™
Dimensions (W x D x H)	260 x 435 x 500 mm
Weight	23 kg
Power consumption	60 W
Frequency	50 / 60 Hz
Connection voltage	100 - 240 VAC ± 10 %
Max. power for all USB-Ports	5 W
IP Code	IP69
Overvoltage category	II
Pollution degree	2
Appliance classes	I
Detector NIR	Thermoelectrically cooled InGaAs
Detector VIS	Si
Wavelength range NIR	900 - 1700 nm
Resolution NIR	7.0 nm
NIR Data Resolution	3.1 nm
Wavelength range VIS	400 - 900 nm
Resolution VIS	Better than 15 nm
VIS Data Resolution	2 nm
Up view illumination spot size	8 mm
Down view illumination spot size	30 mm
Approval	CE / CSA
Lamp type	Tungsten-halogen
Average life (lamp)	7000 h
Display	10.4 in
Storage temperature	max. 45 °C

Ambient conditions

For indoor use only.

Specification	Value
Max. altitude above sea level	2,000 m
Ambient and storage temperature	5–40 °C
Max. relative humidity	80% for temperatures up to 31 °C decreasing linearly to 50% relative humidity at 40 °C

Materials

Component	Material of construction
Housing	Steel 1.4301
Seals housing	EPDM 50
Housing	Aluminum with coating
Seals housing	EPDM 70
Glass up view	Sapphire Al ₂ O ₃
Glass down view	Borofloat
Seals glasses	EPDM A 75
Cooler	Aluminum with coating EPOFLON 526/4562
Sealing frame cooler	HD-PE
Seals cooler	EPDM
Display	Glass
Frame display	Steel 1.4301
Seal Display	1K MS-Polymer Köräpop 225
Seals wire module	FKM
Housing USB/Wi-Fi	Steel 1.4301
Cover USB/Wi-Fi	PE-HD
Seals USB/Wi-Fi	Silicon

Description of function

ProxiMate™ is a NIR spectrometer that can be used to determine the concentration of different parameters in food and feed samples in a nondestructive way.

ProxiMate™ is supplied in different versions. Dependent on the version specified ProxiMate™ is either an NIR or combined NIR and visible spectrometer.

The instrument generates a beam of NIR (and visible) light which is focused onto the sample under investigation. Light reflected from the sample is collected and spatially separated by a diffraction element. The diffracted light is directed onto a diode array detector. Signals from the detector are processed and a reflectance spectrum is constructed. This spectrum undergoes further processing to calculate the constituents required.

Data processing

The NIR light interacts with the sample material in different ways, leaving a characteristic fingerprint on the spectrum. Spectra from both liquids and solids can be measured with ProxiMate™. The spectra of solid samples are collected directly, liquid samples require the use of a transreflectance adapter.

Sample presentation options

The choice of sample presentation of ProxiMate™ is optimised for the type of sample under analysis and for the working environment where the instrument is used.

ProxiMate™ can be configured with a choice of sample presentation options: up view and down view configurations.



NOTE

You can capitalise on the advantages of the up view and down view options in a single instrument.

Up view option

The up view option directs and collects light from the underside of the sample. The NIR light passes through the base of a glass petri dish before interacting with the sample under evaluation. Up view measurement has the advantage that a more consistent surface is presented to the ProxiMate™ ensuring accurate measurement output. Glass petri dishes are recommended to enable best performance. Additionally, when used in conjunction with a transreflectance adapter, it is also possible to measure liquids using the up view option.



NOTE

Choose the up view option for the most consistent measurement performance or for measurements of liquids.

Down view option

With down view option light is focussed onto and collected from the top surface of the sample. In areas where glass is prohibited (such as some food production areas), the down view mode offers the advantage that NIR light does not interact with the sample container. As plastics have their own NIR spectra, changes in type of dish can influence the measurement output leading to perceived measurement shifts. Use of the down view mode prevents this from occurring. Additionally ProxiMate™ down view also allows the use of large volume sample dishes. These are particularly useful with samples that are inhomogeneous, as the measurement is averaged over a much larger sample area.



NOTE

Choose the down view option for areas where glass is prohibited or where larger sample volumes are required.

Spare parts and accessories

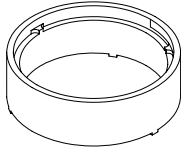
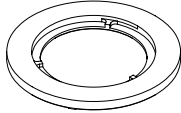
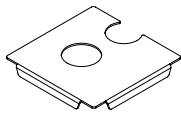
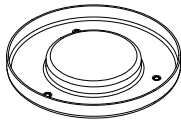
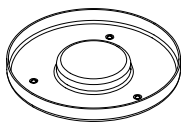
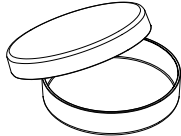
Use only genuine BUCHI consumables and spare parts in order to ensure correct, safe and reliable operation of the system.

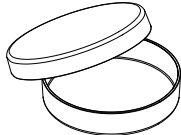
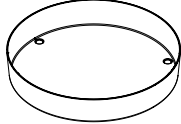
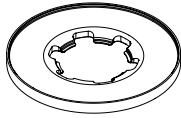
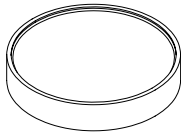
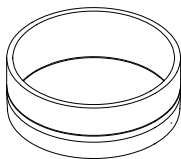
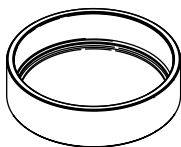
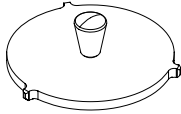
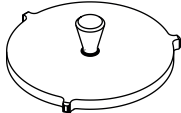
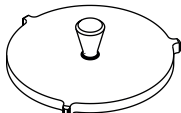


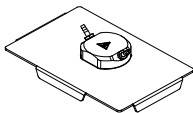
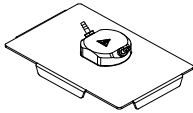
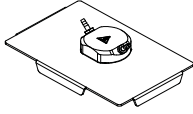
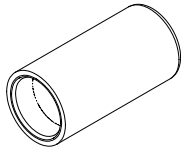
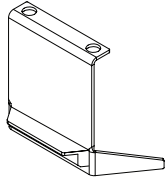
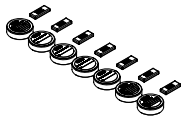
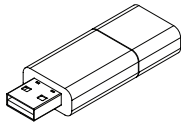
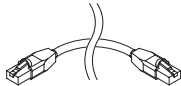
NOTE

Any modifications of spare parts or assemblies are only allowed with the prior written permission of BUCHI.

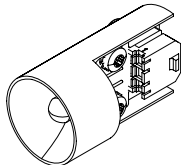
Accessories

	Order no.	Image
External White Reference	11067547	
White reference spacer for down view	11067378	
White reference spacer for up view	11067377	
White reference locating plate	11067391	
PET Large Sample Cup	11067692	
FDA food approved large sample cup	11068015	
Glass petri dishes 10 pcs. (up view) Not suitable for use with Transflectance Cover	11072073	

	Order no.	Image
Plastic petri dishes 240 pcs. (down view)	11066381	
Stainless steel petri-dish For down view measurements only	11074314	
Large sample carrier	11067691	
Carrier for small sample cup (down view)	11065472	
High Performance Sample Cup	11067399	
Robust cup	11055058	
Transflectance cover 0.3 mm Not suitable for use with robust cup	041636	
Transflectance cover 2.0 mm For measurement of crude palm oil. Not suitable for use with robust cup.	11067919	
Transflectance cover for robust cup	11055998	

	Order no.	Image
FlowCell, 0.3 mm For measuring liquid samples under flow conditions	11084993	
FlowCell, 1 mm For measuring liquid samples under flow conditions	11084990	
FlowCell, 2 mm For measuring liquid samples under flow conditions	11084992	
Protection cover USB-WiFi stick	11066582	
Light shield (down view)	11067281	
Performance test standards kit (7 pcs.)	11067545	
NIRWise PLUS Chemometrics suite	11068025	
Set Network cable RJ45 5m	11068780	

Spare parts

	Order no.	Image
Spare lamp	11065441	
Replacement Window HPSC	046246	

Maintenance kits

	Order no.
Kit Customer ProxiMate	11062670
Frequently used wear and spare parts for self-servicing your device, conveniently gathered in one kit	



11594010N en ProxiMate™ 04.2026
Technical data are subject to change without notice
Quality Systems ISO 9001

We are represented by more than 100 distribution partners worldwide.
Find your local representative at:

www.buchi.com

Quality in your hands
